SPECTROSCOPY

ATUL SINGH ARORA



Chemistry Lab III

Dr. K. S. Viswanathan and Dr. K. R. Shamasundar Indian Institue of Science Education and Research, Mohali

August-December, 2012

Every honest researcher I know admits he's just a professional amateur. He's doing whatever he's doing for the first time. That makes him an amateur. He has sense enough to know that he's going to have a lot of trouble, so that makes him a professional.

— Charles F. Kettering (1876-1958) (Holder of 186 patents)

ACKNOWLEDGEMENTS

I express my sincere gratitude to our instructors, Dr. K. S. Viswanathan and Dr. K. R. Shamasundar, for bringing the subject to life and helping us discover, in depth, the science behind the procedures.

I also thank Vivek Sagar (MS11017) for his contribution to this report as my lab-partner, who made the task of performing experiments immensely comfortable and productive at the same time.

CONTENTS

EXPERIMENTS INTRODUCTION TO INFRA RED AND ULTRAVIOLET - VIS-IBLE SPECTROSCOPY 3 1.1 Objective 1.2 Theory 1.2.1 **Basic Concept** 1.2.2 Infra Red Spectroscopy Ultraviolet-Visible Spectroscopy 5 1.3 Experimental Details 1.3.1 Infrared 1.3.2 UV-Vis 1.4 Acknowledgements 2 EXPERIMENT 1A: UV-VISIBLE SPECTROSCOPY 2.1 Objective 2.2 Theory 2.3 Procedure 2.4 Calculations and Measurements 2.5 Spectroscopic Observations and Analysis 9 2.6 Discussion 2.7 Acknowledgements EXPERIMENT 1B: INFRA RED SPECTROSCOPY 13 3.1 Objective 13 3.2 Theory 13 3.3 Procedure 13 3.4 Analysis 13 3.5 Discussion 14 3.6 Acknowledgements EXPERIMENT 2A: JOB'S METHOD OF CONTINUOUS VARI-ATION (UV-VIS) 17 4.1 Objective 17 4.2 Theory 17 4.3 Procedure 17 4.4 Observations and Analysis 18 4.5 Discussion 4.6 Acknowledgements EXPERIMENT 2B: INFRA RED SPECTROSCOPY 21 5.1 Objective 21 5.2 Theory 21 5.3 Procedure 22 5.4 Analysis 22 Reference 5.5 23

 \mathbf{v}

```
5.6 Acknowledgements
6 QUEST FOR SYMMETRY: IR VS. RAMAN SPECTROSCOPY
                                                         25
   6.1 Objective
                   25
   6.2 Theory
                  25
       6.2.1
             Raman Spectroscopy
                                    25
   6.3 Procedure
                    27
   6.4 Observations
                       27
   6.5 Conclusions
                      29
   6.6 Acknowledgements
                             30
   6.7 References
                     30
  ALCOHOL ESTIMATION
                             31
   7.1 Objective
   7.2 Theory
                  31
             Idea
       7.2.1
                     31
       7.2.2
             Details
                        31
   7.3 Procedure
                    32
   7.4 Observations
                       32
   7.5 Acknowledgements
                             32
   7.6 References
                     33
  HYDROGEN BONDING IN ALCOHOL
                                        35
   8.1 Objective
                   35
   8.2 Theory
                  35
   8.3 Procedure
                    35
   8.4 Observations/Discussion
                                  36
   8.5 Acknowledgements
   8.6 References
                     36
  NMR SPECTROSCOPY
                          37
   9.1 Objective
                   37
   9.2 Inspiration
                     37
       Theory
   9.3
                  37
             Chemical Shift
       9.3.1
                               38
             Splitting of Peaks, Integrals
       9.3.2
                                          38
   9.4 Experimental Details
             Locking and Shimming
       9.4.1
                                       39
             Solvent
       9.4.2
             The NMR Jargon
       9.4.3
                                 39
   9.5 Observations/Discussion
                                  39
   9.6 Acknowledgements
   9.7 References
                     40
10 FLOURESCENCE
                     41
   10.1 Objective
                   41
   10.2 Theory
                 41
   10.3 Observations/Discussion
                                  42
   10.4 Acknowledgements
   10.5 References
                     42
11 X-RAY DIFFRACTION AND POLARIMETRY
                                             45
```

11.1 Objective 45	
11.2 Theory 45	
11.2.1 X-ray diffraction	45
11.2.2 Polarimetry 46	
11.3 Observations/Discussion	47
11.3.1 X-ray diffraction	47
11.3.2 Polarimetry 47	
11.4 Acknowledgements 47	
11.5 References 47	
II THE SHOWCASE 49	
BIBLIOGRAPHY 51	

LIST OF FIGURES

T.	NC 1 1 T . C
Figure 1	Michelson Interferometer 4
Figure 2	Least Square Fit 11
Figure 3	Absorbance with Varying Concentrations 20
Figure 4	ATR principle 22
Figure 5	Oxidation of Alcohol 31
Figure 6	Dichromate's Absorption as a function of Al-
	cohol Concentration 34
Figure 7	Hydrogen Bonds in Methanol 36
Figure 8	Concentration vs. Intensity 43
Figure 9	Optical Rotation 48

LIST OF TABLES

Table 1	For the Stock Solution of Benzoic Acid 8
Table 2	For varying concentrations of Benzoic Acid So-
	lution 8
Table 3	Wavelengths Absorbed 9
Table 4	Absorption for various Concentrations 9
Table 5	IR peaks for Urea 14
Table 6	IR peaks for Benzoic Acid 14
Table 7	IR peaks for Thio-Urea 14
Table 8	IR peaks for Naphthalene 15
Table 9	IR peaks for alpha Naphthalene 15
Table 10	Concentrations for Job's Method 18
Table 11	Absorbance for Varying Concentrations 19
Table 12	IR peaks for a CD Cover 22
Table 13	IR peaks for Teflon 23
Table 14	IR peaks for Benzaldehyde 23
Table 15	IR peaks for Phenol 23
Table 16	IR peaks for Benzoic Acid 24
Table 17	IR peaks for Urea 24
Table 18	Pthalic Acid Raman Spectra 28
Table 19	Pthalic Acid IR Spectra 28
Table 20	Pthalic Anydride Raman Spectra 28
Table 21	Pthalic Anydride IR Spectra 29
Table 22	Napthalene Raman Spectra 29
Table 23	Napthalene IR Spectra 29

Table 24	Concentration vs. Intensity	42
·	-	
LISTINGS		

ACRONYMS

Part I EXPERIMENTS

1

INTRODUCTION TO INFRA RED AND ULTRAVIOLET - VISIBLE SPECTROSCOPY

August 9, 2012

1.1 OBJECTIVE

To

- 1. understand the method of Fourier Transformed Infrared and Ultraviolet spectroscopy
- 2. understand the preparation substances for analysis
- 3. handle and operate the Spectroscopes provided

1.2 THEORY

1.2.1 Basic Concept

To find the presence of elements and/or compounds within a given substance, we can use spectroscopy techniques, specially when their concentrations are small and they satisfy certain requirements. The essential idea behind this measurement comes from the fact that elements/compounds absorb lights of certain frequencies to get to a higher energy state. These frequencies are mostly discrete as they correspond to quantized energy levels. This energy could be absorbed for, say, changing the vibrational energy (IR Spectroscopy) or for exciting an electron in the substance to a higher energy level (UV-vis Spectroscopy). We note here that these quantized energy levels are properties of individual substances and are, for most practical purposes, unique.

For the analysis to be possible, the first condition is that the substance must *absorb* light incident to it. Granted this, we can obtain an absorption spectrum for the given substance, which behaves like a fingerprint of the substance. This can thus be used to not only identify the compound, but also to quantify it. For identification, in the simplest case, we simply need to observe the frequency corresponding to the peaks in the absorption spectrum and match it with the known/expected substance(s). Quantification harnesses a rather "obvious" law, termed *Beer-Lambert's Law*. In the simplest form, the law

How much absorption, well, the limit comes from the sensitivity of the experimental setup and concentration of substance given.

quantifies the intuitive notion; higher the concentration of the analyte, higher is the absorption. The relation is given as

$$T = \frac{I}{I_0} = 10^{-\alpha l} = 10^{-\epsilon lc} \tag{1}$$

where I is intensity of incident light, I_0 is intensity of transmitted light, ε is molar absorbtivity, I is the optical path length, and c is molar concentration.

1.2.2 Infra Red Spectroscopy

Infra Red spectroscopy usually deals with energies of the level that cause change in vibrational energies. The wavelength ranges from $2.5\mu m$ to $25\mu m$ ($4000cm^{-1}$ to $400cm^{-1}$). These energies are characteristic for different bonds which is how, using the spectrum, we can identify (and quantify) the bonds present and thus the compound.

The way the spectroscope works for Infra Red, is rather interesting and ingenious. The setup uses a Michelson interferometer.

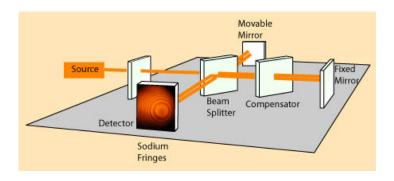


Figure 1: Michelson Interferometer [?]

Before complicating things, let us assume that the source of light is coherent and monochromatic. Now in the Figure 1, we assume the source of light to be the light transmitted through the sample to be analysed. Say the interference at the detector is, at the given configuration, constructive. If we move the movable mirror by $\lambda/4$, where λ is the wavelength of the light, then the detector will receive a dark, destructive interference. If we plot the intensity at the detector as a function of displacement of the moveable mirror, we will, in this case, receive a sine wave.

Now let us crank it up a notch. Let us consider the light to still be coherent, but not monochromatic. Let the source of light contain all the transmitted frequencies. Now if the intensity is plotted against the displacement of the moveable mirror, we will get a superposition of sine waves, and we already know how to decompose them to find individual frequencies using *Fourier Transformation*. This essentially

Coherent
means phase locked
and Monochromatic
means single
frequency

Note the fact that here, the measurement is simultaneous!

gives us the spectrum, with wavenumber (dimensions of one over distance) on one axis, and intensity on the other.

The source of infra-red light is a Silicon Carbide rod, that is heated to produce the desired radiations. The detector is capacitive, or so we were told.

How is SiC coherent?

1.2.3 *Ultraviolet-Visible Spectroscopy*

Here the wavelengths range from 200nm to 900nm. This range is often in the range of differences of electronic energy levels. The light source is a Tungsten-Neon lamp. In this technique, usually movable diffraction grating is used to split the transmitted light into its component frequencies, "one" of which goes to an intensity detector. Which component is detected depends on the orientation of the diffraction grating. Using this information, an intensity vs. frequency plot is generated which is the desired spectrum.

It must be noted here that the process of moving the orientation of the grating, requires time, which means that the measurement for various frequencies' intensities is not simultaneous, which in certain cases might not be suitable. One solution to such a problem is the use of multiple (arrayed) detectors. Yet it doesn't explain why a Michelson interferometer can't be used in this case. The answer to this comes from the fact that the wavelength in this case is much smaller, thereby complicating the task for preparing reflectors.

Wouldn't the ordinary mirrors for "visible light" suffice?

1.3 EXPERIMENTAL DETAILS

For this session, both IR and UV-visible spectroscopy techniques were demonstrated to us in groups of two. For IR spectroscopy, we used KBr as the base and Benzoic Acid as the analyte. For UV-vis, we used water as the base and $KMnO_4$ as the analyte.

1.3.1 Infrared

Procedure

The general procedure is as follows:

- A. Weighed 100mg of KBr and ground it well.
- B. Transferred it the dye (palletizing module).
- c. Converted it into a pellet using the Hydraulic Press (about 8 tons of pressure).
- D. Did IR spectroscopy.
- E. Weighed 100mg of KBr and with it, 2mg of sample (in this case Benzoid Acid) and repeated from item B.

Details about the IR spectroscopy are as follows:

- A. Ran the background to get rid of noise.
- B. Did the accumulation 16 times (set the option).
- c. Subtracted the first spectroscope from the second

Observation

Detected the carboxylic group (peaks for OH and C = O)

1.3.2 UV-Vis

Procedure

The general procedure is as follows:

- A. Took a suitable amount of water.
- в. Did UV-Vis spectroscopy.
- c. Again took water and added a suitable amount of KMnO₄ and mixed it well.
- D. Did UV-Vis spectroscopy.

Details about the UV-vis spectroscopy is corresponding to the previous section.

Observation

Peaks were detected in the range 520nm-550nm, complimentary to Purple as expected.

1.4 ACKNOWLEDGEMENTS

I thank the PhD students, Ms. Shruti and Ms. Shilpa, who helped us with the session and performed the experiment. They also walked us through the procedure for using the apparatus.

August 16, 2012

2.1 OBJECTIVE

To

- 1. prepare 10, 25, 50 and 75 ppm solutions of Benzoic Acid in water
- 2. use one of these for obtaining a spectrogram
- 3. use three of these for calibration and find the concentration of the third experimentally

2.2 THEORY

This section is the same as and has been covered in the previous experiment.

A few simple relations must be stated however, for clarity. Transmissivity or Transmission (T) is related to observables by Equation 2, where I is intensity of the light transmitted, I_0 is intensity of incident light, ϵ is extinction coefficient, l is the path length and c is concentration of the solution.

Relation between Absorbance (A) and intensity of light is as given by Equation 3. Clearly, A has a linear relation with ϵ as given by Equation 4 and we'll use this very relation to find the concentration of the 'unknown' sample.


$$T = \frac{I}{I_0} = 10^{-\epsilon lc} \tag{2}$$

$$A = -\log 10(\frac{I}{I_0}) \tag{3}$$

$$A = \epsilon lc \tag{4}$$

2.3 PROCEDURE

1. First calculated the mass of Benzoic Acid required for making a 25mL a 250 ppm solution in water, to be 6.25 mg.

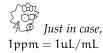

Mass of salt	0.0066	g
Mass of paper	0.0001	g
Mass of Benzoic Acid	6.5	mg
Concentration of Benzoid Acid Used	$\frac{6500}{25} = 260$	ppm

Table 1: For the Stock Solution of Benzoic Acid

10 ppm	$\frac{50}{260} = 0.192$	
25 ppm	$\frac{125}{260} = 0.480$	mL
10 ppm	$\frac{250}{260} = 0.961$	mL
10 ppm	$\frac{50}{375} = 1.442$	mL

Table 2: For varying concentrations of Benzoic Acid Solution

- 2. Weighed roughly the calculated amount of Benzoic Acid and transferred it to a volumetric flask of 25mL, and made the volume 25mL. This solution is henceforth referred to as the stock solution.
- 3. Calculated the volume of the stock solution required to prepare 10, 25, 50 and 75 ppm, 5 mL Benzoic Acid solutions.
- 4. Using a graduated pipette, measured the precise volume in accordance with the calculations and transferred it to a 5mL volumetric flask, and filled it to 5mL of water. Repeated this step for all four concentrations and labelled the flasks accordingly.
- 5. Calibrated the spectroscope with water first to suppress noise, if any.
- 6. Used one of the concentrations to generate a spectrogram.
- 7. Identified the maximum wavelength for use in Beer-Lambert law.
- 8. Used three of the concentrations to calibrate and fourth as an unknown, using the spectroscope.
- 9. Printed the spectrogram and noted the observations off of the screen for the known and unknown concentrations.

2.4 CALCULATIONS AND MEASUREMENTS

Refer to Table 1 for the stock solution and to Table 2 for the variants.

Used the upper meniscus consistently for volume measurements

But does that mean it wouldn't affect the result?

	wavelength (λ)	ABSORBANCE
544.00nm		1.710
524.00nm		1.915
	Table 3: Wavelengths Absorbed	
	CONCENTRATION (c)	ABSORBANCE
75ppm		1.902
75ppm 50ppm 10ppm		0.971
10mm		

Table 4: Absorption for various Concentrations


2.5 SPECTROSCOPIC OBSERVATIONS AND ANALYSIS

The maximum wavelength $(\lambda_{m\alpha x})$ was found to be 524.00nm. Refer to Table 3. The Absorbance (A) for various concentration is given in Table 4

Absorbance of 'unknown' sample (25 ppm): 0.489

In accordance with the application, the best fit for the three points was not linear and the concentration of the unknown sample was calculated to be 30.0720 ppm which amounts to a 20% error.

In accordance with least square straight line fit, we get the slope as 0.0249395 \pm 20.34 % but the intercept as -0.0856124 \pm 310.2%. In accordance with these numbers, the concentration works out to be 16.174 ppm, which is 35% off. Refer to Figure 2 for details.

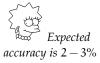

2.6 DISCUSSION

The accuracy of this experiment, as told by our instructor, is expected to be less than 3%. Since we clearly failed at achieving that number, we analyse in this section, what could've caused the enormous error.

The error could be caused by making errors in any of the following processes

- 1. Weighing
- 2. Volume measurement
 - a) Pipetting
 - b) Filling the Volumetric Flask

Despite taking sufficient care, even if one of the steps mentioned above go wrong, for even one of the samples, since the number of samples for calibration is just three, the entire experiment can, and in this case, did go wrong.

Thus for ensuring such errors do *not* crop up in the future, we can attempt the following

1. Weighing:

- a) Put the paper on the weighing scale and tare it
- b) Add the required amount of substance ensuring nothing gets wet
- c) Keep the container, viz. Volumetric Flask, with accessories, viz. Funnel etc. near the setup
- d) After transferring the contents, weigh the paper again to ensure everything got transferred

2. Volume Measurement


- a) Pipetting
 - i. Read the upper meniscus consistently (for this experiment and in general for opaque solutions)
 - ii. Use graduated pipettes that can provide sufficient precision
 - iii. Ensure the liquid doesn't stick on the walls
- b) Filling Volumetric Flask
 - i. Again, use the upper meniscus
 - ii. Use a dropper to add towards the end to avoid accidentally adding extra volume

2.7 ACKNOWLEDGEMENTS

I thank Ms. Sumyra, the PhD student, who helped us with the use of the spectroscope. I also acknowledge the contribution of Ms. Athira J. Nair and Mr. Arpit Porwal for the performance of the experiment as team members.

This is most likely what went wrong!

Experiment: UV-visible Spectroscopy

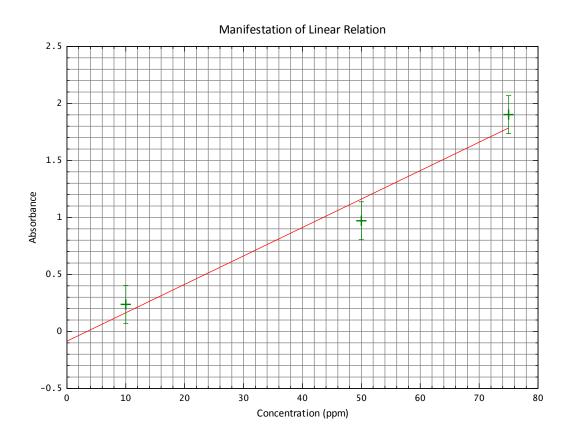

Slope of Best Fit Line : +0.0249 Intercept of Best Fit Line : -0.0856

Performed on: August 16, 2012 Performed by: Athira, Arpit and Atul

Figure 2: Least Square Fit

EXPERIMENT 1B: INFRA RED SPECTROSCOPY

August 23, 2012

3.1 OBJECTIVE

To study the spectroscopic characteristics of atleast three of the following compounds using an FTIR spectroscope.

- 1. Benzoic Acid
- 2. Naphthalene
- 3. Urea
- 4. Thio-Urea
- 5. Salicylic Acid

3.2 THEORY

This section is the same as and has been covered earlier in Chapter 1.

3.3 PROCEDURE

- 1. Weighed approximately 200mg of KBr and 2-4mg of sample and ground it well using mortar and pestle. (Didn't use any sample for the first time, to use for baseline correction.)
- 2. Transferred the powdered contents into the palette maker and put it inside the hydraulic press
- 3. Applied about 6-8 tons of pressure for 30-60 seconds. Released the pressure after, and transferred the contents to the palette holder.
- 4. Placed the palette holder into the FTIR spectroscope and scanned it (in case of the first palette, background scanned it).

3.4 ANALYSIS

- 1. Urea Table 5
- 2. Benzoic Acid Table 6

The process must be done quickly to avoid accumulation of moisture.

N – H stretching	3399.69	cm^{-1}
C = O stretching	1669.64	cm^{-1}
N – H bending	1625.39	cm^{-1}
	Table 5: IR peaks for Urea	
C – H stretching	2832.14	cm^{-1}
O – H stretching	2554.86	cm^{-1}
C = O stretching	1694.66 and 1699.26	cm^{-1}

Table 6: IR peaks for Benzoic Acid

- 3. Thio Urea Table 7
- 4. Naphthalene Since there are no functional groups, according to L. G. Wade's *Organic Chemistry*, peaks are expected to be observed in the following ranges
 - a) For Naphthalene in general:
 - i. between 3000 and 3100
 - ii. between 1570 and 1650

Corresponding observations are given in Table 8

- b) For alpha naphthalene, we have
 - i. between 1375 and 1425
 - ii. between 750 and 810

Corresponding observations are given in Table 9

- c) For beta naphthalene, peaks are expected to be in the ranges:
 - i. less than 700
 - ii. between 800 and 860

These were found to be absent.

3.5 DISCUSSION

For this experiment, even though we knew which compounds we're taking, after printing the spectrographs, we didn't label them to see if

NH ₂ stretching	3368.14	cm^{-1}
NH ₂ bending	1588.41	cm^{-1}
C = S stretching	1092.87	cm^{-1}

Table 7: IR peaks for Thio-Urea



3061.97		cm ⁻¹
3048.41		cm^{-1}
3029.35		cm^{-1}
1651.94		cm^{-1}
1634.25		cm^{-1}
1592.53		cm^{-1}
	Table 8: IR peaks for Naphthalene	
1383.98		cm ⁻¹
782.47		cm ⁻¹

Table 9: IR peaks for alpha Naphthalene

they could be identified by matching them with their expected peaks. This method does in fact work quite accurately and the fingerprint region was also observed to match, when compared with known spectrographs. I have still not been able to understand the meaning of alpha and beta structures of Naphthalene, which I so far believe to not exist. However, experimentally, in accordance with the literature, the structure seems to be alpha.

3.6 ACKNOWLEDGEMENTS

I acknowledge the contribution of Mr. Arpit Porwal for the performance of the experiment as a team member. I would also like to thank Mr. Biplob Nandy and Ms. Saumya Gupta for exchanging results to mutually provide scope for wider analysis.

I am grateful to Mr. Arjit Kant Gupta who helped me with the analysis and confirmation of the results by sharing his spectrographs of known compounds.

I also thank L. G. Wade for the book "Organic Chemistry" which was used to look up the expected IR peaks in the compounds analysed.

EXPERIMENT 2A: JOB'S METHOD OF CONTINUOUS VARIATION (UV-VIS)

August 30, 2012

4.1 OBJECTIVE

To find the reaction stoichiometry of the Fe^{3+} —Salicylic complex using Job's method of continuous variation.

4.2 THEORY

The idea is fairly simple. We analyse a compound which is formed by a combination of two reactant substances. The only requirement is that we should know of some way using which we can quantify the compound, given a mixture consisting of the compound along with the residual reactant substances.

Here's what we do...:

- 1. Decide on a total number of moles you will initially take (of both the reactant substances combined)
- 2. Now take the reactant substances in various ratios, such that the total number of moles is as decided
- 3. For each ratio, find out the quantity of the compound obtained ...and why:

A little thought will make the entire procedure appear elegantly simple. In the various combinations, one of the reactants will always be a limiting reagent. However, the highest yield will be obtained in the case when the reactants are closest to the stoichiometric ratio of the compound. Thus, a plot of the concentration of the compound against that of one of the reactants will attain a maxima. The point at which the maxima is obtained, can be used to find the concentration of both reactants, whose ratio gives us the reaction stoichiometry, with respect to the reactants!

And just to state the obvious, as the title suggests, the technique we use for quantifying is UV Visible spectroscopy.

4.3 PROCEDURE

As always, when practically performing an experiment, we have a few extra details to take care of.

Fe^{3+} solution's volume (mL)	SALICYLIC ACID'S VOLUME (mL)
0.5	4.5
2.0	3.0
2.5	2.5
3.0	2.0
4.5	0.5

Table 10: Concentrations for Job's Method

- 1. Make a 0.001M, 25mL solution of Iron Nitrate, $Fe(NO_3)_3$.
 - a) The Molar Mass was given to be 404g.
 - b) Thus, for 25mL we need 10.0mg.
- 2. Make a 0.001M, 25mL solution of Salicylic Acid
 - a) The Molar Mass was given to be 138.12g
 - b) Thus, for 25mL, we need 3.5mg.
- 3. We use 10mL as the total volume.
- 4. We used the volumes given in Table 10, of the reactants from the 0.001M solutions prepared, for creating 5mL solutions and marked the volumetric flasks with the corresponding concentrations.
 - a) Used an appropriate graduated pipette for measuring the volumes
 - b) Used the pipette for measuring one of the reactants only, and filled the volumetric flask with the other reactant using a dropper, to 5mL using the mark on the flask. This was done to avoid errors.
 - c) Each of the flasks were ensured to be dry.
- 5. Recorded the spectrum of one of the samples and found the peak corresponding to the compound we're interested in, viz. Fe^{3+} —Salicylate Complex.
- 6. Recorded the absorbance of the frequency determined in the previous step for all concentrations.

4.4 OBSERVATIONS AND ANALYSIS

- 1. Characteristic wavelength was found to be 518.00.
- 2. Intensities for varying concentrations of Fe^{3+} have been listed in Table 11.

Since the molarity of both solutions, is the same, we can use volume as a measure of number of moles.

	Fe^{3+} solution's volume (mL)	ABSORBANCE
0.5		0.350
2.0		0.682
2.5		0.801
3.0		0.740
4.5		0.282

Table 11: Absorbance for Varying Concentrations

3. From Figure 3, it's clear that the maximum concentration is obtained for the *stoichiometric ratio of 1:1*, as the maxima is obtained very close to 2.5mL, which represents equal concentrations of both reactants.

4.5 DISCUSSION

Fortunately or unfortunately, the experiment returned expected results and therefore there's no real requirement of a discussion. However, it must be mentioned that by just looking at the colour of the solution, the search can be narrowed to a great extend, that is, cases in which the product and/or reactants are coloured.


4.6 ACKNOWLEDGEMENTS

I acknowledge the contribution of Mr. Arpit Porwal and Ms. Athira Nair, for the performance of the experiment as team members. I would also like to thank Ms. Ritu Roy Chowdhury for discussing the expected graph, which helped me understand the experiment better.

How do you know the graph will be linear on both sides?

[Product] α [Limiting Reactant]

Experiment: Job's Method

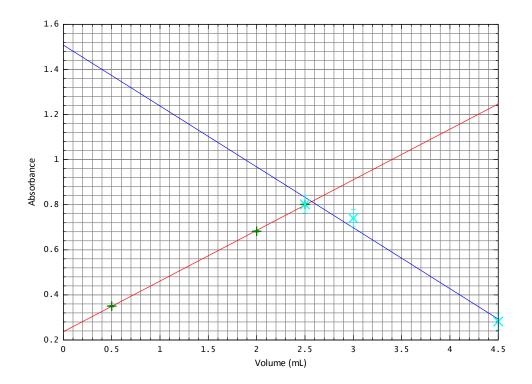

Slope of Best Fit Line 1 : +0.2245 Standard Deviation 1 : +0.0035 Slope of Best Fit Line 2 : -0.2705 Standard Deviation 2 : +0.0385

Performed on: September 3, 2012 Performed by: Athira, Arpit and Atul

Figure 3

September 13, 2012

5.1 OBJECTIVE

To study the spectroscopic characteristics of the following compounds using an ATR-FTIR spectroscope.

- 1. Polymers
 - a) CD Cover
 - b) Teflon
- 2. Solids
 - a) Benzaldehyde
 - b) Phenol
- 3. Liquids
 - a) Benzoic Acid
 - b) Urea

5.2 THEORY

This FTIR part of this section is the same as and has been covered earlier in Chapter 1. What makes this section special is the word ATR, which stands for Attenuated Total (internal) Reflection. Total internal reflection occurs, as we know, when light travelling in an optically denser medium, is incident on a rarer surface, with an angle greater than some 'critical' value, then the light is completely reflected. What has this got to do with spectroscopy, well as it turns out, light when is reflected as aforesaid, it penetrates to about 20 micrometers into the rarer surface. Here's where the fun begins. As the wavelength approaches the absorption range of the rarer substance, the reflected ray undergoes higher attenuation. And that does it! We harness this very property to 'scan' through different wavelengths and find the spectrum of the substance.

Experimentally, we use a crystal of a substance of very high refractive index, typically Thallium Bromide/Iodide and let the incident light from the FTIR, enter the crystal as shown in Figure 4, hit the rarer surface multiple times to accumulate attenuation and then allow this light ray to pass through the same FTIR detector. (In our case however, we used ZnSe.)

C — H stretching	2917.68	cm^{-1}
C — H bending	1376.45	cm^{-1}

Table 12: IR peaks for a CD Cover

What is the immediate advantage? Aside from the ability to measure spectra of opaque substances, we can use this technique to find spectra of fluids as well. Further, since this is a surface technique, situations where bulk is not of interest, this technique comes to the rescue.

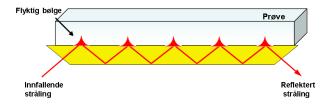


Figure 4

5.3 PROCEDURE

The procedure for ATR is rather straightforward and easily done, unlike the palette preparation requirement in the conventional FTIR.

- 1. Method for Liquids:
 - a) Cleaned the ZnSe tray with CCl₄ and wiped off its excess
 - b) Liquid samples were then loaded and the tray covered
- 2. Method for Solids:
 - a) Powdered the solid finely using Pestle and Mortar and spread evenly on the tray.
 - b) The tray was placed and tightened using a screw firmly.
- 3. For Films and Plastics:
 - a) An extra apparatus provided with the ATR crystal was used.

5.4 ANALYSIS

- 1. CD Cover Table 12
- 2. Teflon Table 13
- 3. Benzaldehyde Table 14

C — F stretching	1206.97 and 1150.19	cm^{-1}
Table 13: IR peaks for Teflon		
O — H due to presence of Benzoic Ac	id 3059.5	cm^{-1}
H - C = O stretch	2738.04	cm^{-1}
Carbonyl $C = O$ stretch	1583.77	cm^{-1}

Table 14: IR peaks for Benzaldehyde

- 4. Phenol Table 15
- 5. Benzoic Acid Table 16
- 6. Urea Table 17

5.5 REFERENCE

- 1. Modern Spectroscopy, Fourth Edition, by J. Michael Hollas
- 2. http://en.wikipedia.org/wiki/Total_internal_reflection
- 3. Lecture Notes from CHM211 and CHM201, 2012-13
- 4. Organic Chemistry by L. G. Wade

5.6 ACKNOWLEDGEMENTS

I acknowledge the contribution of Mr. Arpit Porwal and Ms. Athira Nair, for the performance of the experiment as team members. I specially thank Mr. Arpit Porwal for assisting me with the print outs.

C – H stretch	3049.49	cm^{-1}
O – H stretch	3324.03	cm^{-1}
C – C multiple bond stretching	1450 - 1600	cm^{-1}
C – O stretching vibration	1367.17	cm^{-1}
O – H bending	1222.87	cm^{-1}
C – H bending	748.40	cm^{-1}

Table 15: IR peaks for Phenol

C – H stretch	2822.7	cm^{-1}
O – H stretch	2554.7	cm^{-1}
C = O stretch	1688	cm^{-1}

Table 16: IR peaks for Benzoic Acid

C = O stretch	1598.07	cm^{-1}
C - N stretch	1460	cm^{-1}

Table 17: IR peaks for Urea

QUEST FOR SYMMETRY: IR VS. RAMAN SPECTROSCOPY

September 20, 2012

6.1 OBJECTIVE

To gain an insight into the structure of the following chemicals by analysing their Raman and IR spectra:

- 1. Pthalic Acid
- 2. Pthalic Anhydrate
- 3. Napthalene

6.2 THEORY

6.2.1 Raman Spectroscopy

Idea Behind it

What's so cool about Raman Spectroscopy? Well the fact that it uses scattering of light instead of absorption. This has a remarkable consequence. Let's start from the beginning. There are various numbers of molecules in any given quantum state, frequency of which depends on Temperature amongst various other factors (ambiguity is necessary for generality). Now, when an incident beam, of a given frequency, say v_i , strikes a molecule, it is excited from whichever allowed quantum state it was in, to an energy state, called a *virtual state*, (which in this case is higher in energy). Now the speciality of a virtual state is that it can exist anywhere in the energy scale, however its lifetime tends to zero. They are not 'physical', whatever that's supposed to mean! So immediately after jumping to the virtual state, the molecule drops to an allowed energy state, where by allowed, I mean the states allowed by the selection rules.

Stokes and Anti-Stokes

Now say the energy difference between the two energy states, viz. final minus initial, is represented by $\nu_{absorbed}$. Energy of the scattered photon is given by

$$v_{s} = v_{i} - v_{absorbed} \tag{5}$$

It's then obvious to note that if $\nu_{absorbed}$ is positive, then the scattered photon will have a lower energy, and when $\nu_{absorbed}$ is negative, the energy is higher. But here's something which is not immediately obvious. When the resultant energy is lower, corresponding spectroscopic lines are called stokes lines and in vibrational-raman spectroscopy, they are brighter than their counterpart, the anti-stokes lines (definition of which is implied). However in case of rotational-raman, the intensities of both the anti-stokes and stokes lines are comparable.

Vibrational Raman Spectra

Let us lose the generality and talk about, Vibrational-Raman Spectroscopy. We can borrow the expression for energy levels (expressed in wavenumbers) from the prior discussions as:

$$\varepsilon = \overline{\omega}(\nu + \frac{1}{2}) - \overline{\omega}_e x_e (\nu + \frac{1}{2})^2 \text{ cm}^{-1}$$
 (6)

where the symbols have the usual meanings. The selection rule for this reads

$$\Delta v = 0, \pm 1, \pm 2, \dots \tag{7}$$

Transitions from $v=0 \Rightarrow v=1$ are fundamental and are clearly visible on the spectra due to high population of molecules in the state defined by v=0. Transitions from $v=1 \Rightarrow v=2$ (Hot Bands) and $v=0 \Rightarrow v=2$ (overtones) are not observable in most cases. This can be attributed to low populations. Thus, since $v_{absorbed}=v_{fundamental}$, we can expect a spectra determined by Equation 5, with Stoke's lines showing up with a frequency less than v_i and Anti-Stokes' mirrored about v_i with a much lower intensity.

What does Raman 'see'?

As Raman spectroscopy is essentially based on scattering of light, whether the scattered light will be scattered 'elastically' or 'inelastically' would depend on whether change in molecular configuration results in a change in polarisability. This is best understood by considering a two simple examples.

1. Consider a Carbon Dioxide molecule and imagine it executing a symmetric stretch. Now the polarisability of the molecule will increase as the oxygen atoms move away from the equilibrium position. Contrary to this, the polarisability will drop when the oxygen atoms congest with the carbon. This essentially means if we plot polarisability, α against displacement from equilibrium position of the oxygen atom, ζ , the slope at $\zeta=0$ is positive. Note this. Also note, the dipole moment doesn't change in this case.

What about the overtones then?

Note polarisability is distinct from the dipole moment IR is sensitive to

2. Again consider a Carbon Dioxide molecule but this time, imagine its bending, with the Carbon in the centre and Oxygens oscillating up and down. In this case, its easy to observe that the polarisability will be least when the molecule is its equilibrium condition. Thus, in this case, slope of the α vs ζ graph, at $\zeta=0$ will be zero. Also observe, the dipole moment changes as a function of ζ .

Thus, Raman will 'see', in this example, only the first, viz. Symmetric Stretch, which doesn't show up in IR. And as you would've guessed by now, IR will 'see' the second, but Raman won't.

Consequences of Symmetry

Clearly symmetry has a major role to play in deciding what characteristic of the molecule shows up where. Discussing that at length is not in the scope of this report. Yet the following must be stated, for the analysis relies on understanding of the principle.

"Rule of mutual exclusion ¹: If a molecule has a *centre of symmetry*, then Raman active vibrations are infra-red inactive and vice versa. If there is no centre of symmetry then some (but not necessarily all) vibrations may be both Raman and infra-red active."

The converse of the statement is also true, viz. if the molecule shows no common lines in Raman and infra-red, then the molecule has a centre of symmetry (although concluding this might be tricky as the peaks can be diminishingly small).

6.3 PROCEDURE

For ATR, please refer to Chapter 5. For Raman, we performed two calibration follow these simple steps

- 1. The sample was taken in a glass slide
- 2. The glass slide was placed inside the Raman apparatus
- 3. Using the attached computer, the point for analysis was fixed
- 4. Spectrum was recorded

6.4 OBSERVATIONS

- 1. Pthalic Acid For Raman, refer to Table 18 and for IR, refer to Table 19
- 2. Pthalic Anhydride For Raman, refer to Table 20 and for IR, refer to Table 21

Wait, what is ζ supposed to be here?

It's the angle of Carbon Oxygen bond from the equilibrium position

¹ Fundamentals of Molecular Spectroscopy, 4th Edition, Banwell and McCash

PHENOMENON	EXPECTED (CM^{-1})	OBSERVED (CM ⁻¹)
C – H Out of plane	690 — 900	773
bending		
C – O stretching	1000 - 1300	1047, 1177
C = C stretching	1475 - 1600	1644
C — H stretching	3050 - 3150	3092

Table 18: Pthalic Acid Raman Spectra

PHENOMENON	EXPECTED (CM^{-1})	OBSERVED (CM^{-1})
C – H Out of plane	690 - 900	794.16, 739.58
bending		
C - O stretching	1000 - 1300	11071.27, 1279.40
C = C stretching	1475 - 1600	(aprox) 1530
C – H stretching (over-	2400 - 3400	2871.26
laps with O – H)		

Table 19: Pthalic Acid IR Spectra

PHENOMENON	EXPECTED (cm^{-1})	OBSERVED (CM^{-1})
С — Н	690 - 900	773, 737
Out of		
plane		
bend-		
ing		
C - O	1000 - 1300	1047, 1009 (only strong listed)
stretch-		
ing		
C = C	1475 - 1600	1600
stretch-		
ing		
C – H	3050 - 3150	3081
stretch-		
ing		
C = O	1800 - 1830,1740 - 1775	1848, 1766
stretch-		
ing		

Table 20: Pthalic Anydride Raman Spectra

PHENOMENON	EXPECTED (CM^{-1})	OBSERVED (CM^{-1})	
C – H Out of	690 - 900	712	
plane bending			
C - O stretching	1000 - 1300	1282	
C = C stretching	1475 - 1600	(absent)	
C — H stretching	3050 - 3150	3023 (slightly out)	
C = O stretching	1800 - 1830,1740 - 1775	5 1820, 1764	

Table 21: Pthalic Anydride IR Spectra

PHENOMENON	EXPECTED (CM^{-1})	OBSERVED (CM^{-1})
C – H bend	690 — 900	764
C = C stretch	1475 - 1600	1578
C = H stretch	3050 - 3150	3062

Table 22: Napthalene Raman Spectra

3. Napthalene - For Raman, refer to Table 22 and for IR, refer to Table 23

6.5 CONCLUSIONS

- 1. Acid and Anhydride
 - a) The energy levels split in case of Anhydride because the C – OH get very close compared to their position in the Acid. The split can be attributed to in phase and out of phase motion. Further the motions get coupled in case of Anhydride.
 - b) The peaks obtained were not complimentary, therefore we could conclude both molecules lack a centre of symmetry.

2. Naphthalene

a) It should be silent in IR, however we observed peaks. Reason is that the region is not isotropic. The crystal structure distorts the molecule and perturbations are induced, which

PHENOMENON	EXPECTED (CM^{-1})	OBSERVED (CM^{-1})
C – H bend	690 — 900	776
C = C stretch	1475 - 1600	1505, 1590
C = H stretch	3050 - 3150	3049

Table 23: Napthalene IR Spectra

makes the molecule visible in IR, despite having a centre of symmetry.

6.6 ACKNOWLEDGEMENTS

I thank our PhD student guide, Ms. Shruti, who helped us with the performance of the Raman experiments. I thank Ritu Roy Chaudhury, Prashansa Gupta, Athira T John and Vivek Sagar for discussing their analysis.

6.7 REFERENCES

- 1. Modern Spectroscopy, Fourth Edition, J. Michael Hollas
- 2. Fundamentals of Molecular Spectroscopy, 4th Edition, Banwell and McCash
- 3. Spectroscopy, Kris Bybyan

7

September 27, 2012

7.1 OBJECTIVE

To estimate the level of alcohol in a given solution (in water).

7.2 THEORY

7.2.1 Idea

The objective is to measure the concentration of alcohol. So we first think of the methods by which alcohol can be recognized. One such method would be to oxidise the alcohol using a reagent that changes colour. This way we can use the UV-Vis spectrometer to quantify the concentration of the said reagent. One reagent that fits the bill is Potassium Dichromate ($K_2Cr_2O_7$).

7.2.2 Details

We want Potassium Dichromate in its +6 oxidation state. For this, it must have a sufficiently acidic medium for which H_2SO_4 is used. Further, $AgNO_3$ is required for speeding up the reaction. The basic reaction is give in Figure 5.

Figure 5

7.3 PROCEDURE

A 50% (ν/ν) solution of H₂SO₄ was provided as 'Solution A'. A 200 μ L/L solution of Ethanol was provided as 'Solution B'. The following steps were then followed:

- 1. Weighed 20 mg of $K_2Cr_2O_7$ and 25 mg of $AgNO_3$ into a volumetric flask of 50 mL using a funnel. Solution A was used to make the volume to 50 mL. The solution thus obtained was labelled 'Solution C'.
- 2. Pipetted the following volumes in four separate 10 mL volumetric flasks.
 - a) 7.5 mL of solution C and 0.5 mL of solution B
 - b) 7.5 mL of solution C and 1.0 mL of solution B
 - c) 7.5 mL of solution C and 1.5 mL of solution B
 - d) 7.5 mL of solution C and 2.0 mL of solution B

And in another 10 mL flask, pipetted 7.5 mL of solution C and a volume unknown to other team members between 0.5 to 2.0 mL of solution B.

The volume was filled up using Solution A.

3. Measured the absorbance as a function of alcohol concentration and estimated the concentration of unknown alcohol sample given.

7.4 OBSERVATIONS

Our experiment failed since the trend was unacceptable. We will repeat the experiment in the make-up lab, as was advised by our Professor.

November 15, 2012

The experiment was repeated in accordance with the advice and the results are given in Figure 6. We had made a mistake in measurement and addition of solution B in the fourth solution; it hasn't thus been included in the result.

7.5 ACKNOWLEDGEMENTS

I thank our PhD student guide, Ms. Shruti, who helped us with the performance of the UV-Vis Spectroscopy. I also acknowledge the contribution of my team members, Ms. Athira and Mr. Arpit, for performance of the experiment.

November 15, 2012

Even if the mass is not very precise, it doesn't matter so long as the same volume of the solution is used, since the important factor here is change in concentration



I would like to thank Mr. Arpit for performing the experiment again with me. Also, I must acknowledge the contribution of our Lab Assistants who helped us operate the UV-Vis spectroscope. And as always, Prof. KSV is sincerely acknowledged.

7.6 REFERENCES

1. K. S. Visvanathan's Notes

Absorption by Dichromate Ions	Concentration of Ethenol in 200 uL/L
0.5	0.327
1	0.231
1.5	0.15
0.75	0.309

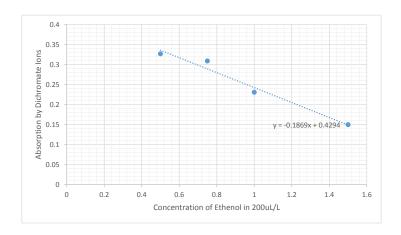


Figure 6

8

HYDROGEN BONDING IN ALCOHOL

September 27, 2012

8.1 OBJECTIVE

To study the hydrogen bonds in alcohols using infrared spectroscopy.

8.2 THEORY

Alcohols in liquid form, interact via hydrogen bonds. This interaction results in a broad peak in the spectrum of alcohols, near the O-H stretching region, viz. 3500 cm⁻¹.

To verify this theory, we dilute the alcohol progressively in CCl₄ and expect the peak to get narrower as fewer interactions take place, with lower concentration. One may argue that the intensity of the peak will also drop. However, this method increases 'contrast' at the expense of 'brightness', viz. the peak relative to the rest of the spectra will continue to be observable and its spread is expected to decrease.

We used Methanol in our experiment. Hydrogen bonds form as shown in Figure 7

8.2.0.1 A Little Deeper

The monomeric alcohol is responsible for producing a sharp peak. The broadening is attributed to the associated alcohol. However, that's not it. Now let's consider a broad peak. Here, the maxima will correspond to the associated molecules, with O-H vibrating with an energy corresponding to that particular wavenumber, where the peak is observed. We can safely conclude that due to association, the molecules will vibrate with lower energy and can therefore expect the peak to be red shifted, with respect to the sharp peak of monomeric alcohol.

Also, its natural to conclude that the broadness and extent of shift signifies the extend of hydrogen bonding and/or association of alcohol under investigation.

8.3 PROCEDURE

We took Methanol (other groups took Ethanol and Propanol because the analysis takes longer, so all groups couldn't afford to 2) and CCl₄ and three 10 mL volumetric flasks. ¹

¹ Reformulated from Prof. KS Viswanathan's Notes

Figure 7

- 1. In the first, pipetted 0.5 mL of alcohol and made up the volume using CCl₄. Labelled this as Solution A.
- 2. In the second, took 5 mL of solution A and made up the volume again using CCl₄. Labelled this as solution B.
- 3. In the third, took 5 mL of solution B and made up the volume as before.



Recorded the IR Spectra of the solutions prepared above, and of neat alcohol . Then compared the spectra for varying concentrations and compared the characteristics of the peak corresponding to hydrogen bonding.

8.4 OBSERVATIONS/DISCUSSION

The graphs and observations have been attached herewith, however they do not seem to be consistent with the thoery. Since K. S. Visvanathan sir was not in the campus until the date of submission, I couldn't understand how to interpret the data.

8.5 ACKNOWLEDGEMENTS

I thank our Lab Assistant, Mr. Mangat Kashyap, who helped us with the performance of the UV-Vis Spectroscopy. I also acknowledge the contribution of my team members, Ms. Athira and Mr. Arpit, for performance of the experiment.

8.6 REFERENCES

1. K. S. Visvanathan's Notes

9

October 11, 2012

9.1 OBJECTIVE

To obtain and analyse an NMR spectrum of the following compounds:

- 1. Benzaldehyde
- 2. Ethyl Acetate

9.2 INSPIRATION

Writing the record for this experiment was a rather unique experience as NMR spectroscopy is amongst the first that we've done practically without having studied the theory, and it doesn't use the conventional methods of absorption of light.

9.3 THEORY

Let's start by talking about the magnetic spins, not of electrons, but of protons instead. We take for granted the fact that the quantization of projection of magnetic spin of a proton is conceptually the same as that of an electron. We note that magnetic moment of the 'spinning' charge, would be proportional to the angular spin. We also observe that there's no way to distinguish between a $+\frac{1}{2}$ and $-\frac{1}{2}$ spin in an isolated nucleus, however if we apply an external constant magnetic field, something interesting happens. The energy of the two spin states, splits into two. To understand how and why this happens, observe that the external magnetic field, will cause the magnetic field of the proton to align along it. However, since the projection of the magnetic field can only be $+\frac{1}{2}$ or $-\frac{1}{2}$ and not zero, thus the magnetic field can't point along the external field, and instead it precesses such that the projection is $+\frac{1}{2}$ along the direction of the external magnetic field. Consequently, the only other possibility is $-\frac{1}{2}$, where the proton's magnetic field vector opposes the external magnetic field. Thus, in accordance with our nomenclature, the $+\frac{1}{2}$ state is lowered in energy, and is thus distinguishable from the $-\frac{1}{2}$ energy state.

The projection is along some fixed arbitrary axis, which here is implicitly assumed to be along the magnetic field for reasons which will become clear soon

¹ This is an intuitive argument, but the end result is consistent with precise calculations

9.3.1 Chemical Shift

As was previously mentioned, the energy difference ΔE in the two states, depends on the Magnetic Field. Now here's where things get further interesting. Protons sit inside the nucleus, with electrons around. Electrons also have their own magnetic field which *shields* the external field from the proton. Thus, assuming we keep the magnetic field constant, the value ΔE *grossly* depends on the electron density around it. Thus we already know a little about the environment of the proton. This phenomenon, viz. change in ΔE for the same nucleus, because of electron distribution, is referred to as *Chemical Shift*. For instance, if we have a CH_3CH_2-OH , then the H of OH will have a higher ΔE compared to the H of CH_3 . This is because oxygen is a highly electronegative atom, which would cause the electron density to be lower in the vicinity of the Hydrogen atom, thereby exposing it to a higher Magnetic field and thus, a higher ΔE .

If they're paired (true for electrons and protons), they're silent. In fact, ∃ something called Electron Spin

Resonance

electrons also talk in

NMR? You've after

all assumed them to

protons, except for the charge of course.

be identical to

9.3.2 *Splitting of Peaks, Integrals*

Now let's us take a finer look at what could affect the magnetic field of a proton. To proceed with this section, let us again consider the example of $CH_3CH_2 - OH$. Consider the proton of the methyl group first. Which proton, well they're indistinguishable. So we would expect to see just one peak corresponding to them (where by peak I mean the peak corresponding to ΔE , how the peak is obtained will be discussed in the next section). However, note that on the adjacent carbon, two more protons are sitting. They can either both be plus half spin, or both be minus half spin, or one plus one minus. Since a proton is roughly equally likely to be either in the plus half state or in a minus half state, we will have molecules in + +, - +, - - proton spin states, in the ratio 1:2:1 respectively (the order of plus minus doesn't matter, viz. + - is the same as - +, which is why the second term is 2 in the ratio). Now here's why this discussion is even worthwhile; the molecules with protons (protons of CH2 are being referred to) in + - spin state, will have zero net magnetic moment, and thus will not affect the ΔE of the protons on the methyl in any way. However, molecules that have protons in + + state will add to the external magnetic field experienced by the protons of the methyl group, thereby increasing their ΔE value. Symmetrically, the - - state will decrease their ΔE value. So, the number of molecules with ΔE unchanged to those with ΔE greater or smaller will be 2:1, in accordance with the ratio obtained earlier. Consequently the peaks would have an area ratio of 1:2:1 as is apparent in the graph. By the same line of arguments, we conclude that the protons in CH₂, will have 1:3:3:1 ² ratio of area

² The concept can be easily generalized to n spin-spin coupling and that gives n+1 as multiplicity, with ratios in accordance with the pascal's triangle.

of peaks, affected by the three protons in the methyl group. Please carefully note that we've assumed the energy level of the neighbouring proton fixed in our analysis and this is a valid approximation since the effect of this split is negligible compared to ΔE of both, the neighbour and affected protons.

9.4 EXPERIMENTAL DETAILS

9.4.1 Locking and Shimming

Since the value of ΔE is highly dependent on the magnetic field, there mustn't be any gradient in the field. To attain this, two processes are used. Shimming reduces the width of the peak, which is achieved by dropping the gradient using various small electromagnets. To keep the position of the peek fixed, a process known as locking is employed, which also, essentially requires stabilization of the magnetic field.

9.4.2 Solvent

The solvent used here is CDCl₃ and it has a deuterium instead of a proton which doesn't talk in the NMR region we're interested in, in this case. Also, the solvent is volatile, which is very helpful since after the experiment, it helps itself out, leaving the sample intact.

9.4.3 The NMR Jargon

Since the value of ΔE depends on the Magnetic Field used, results from different machines become difficult to compare. For this, a simple normalization of sorts is done. From the ΔE of the sample, ΔE of a standard is substracted, which is TMS. It's used for a variety of reasons; one it has only one kind of hydrogen (you don't want your standard to have multiple peaks), two it's peak lies far away from the molecules we usually take interest in, since Silicon is electropositive. So, after subtraction, we divide by the Magnetic Field strength. The number we get comes in the order of 10^6 and this we express as $j \times 10^6$ so that j becomes a part in a million. It is this number that we see on the graph as ppm!

9.5 OBSERVATIONS/DISCUSSION

The spectrum obtained has been appended to this experiment, with the peaks assigned in accordance with the theory and literature, for Ethyl Acetate and Benzaldehyde respectively.

9.6 ACKNOWLEDGEMENTS

I thank the folks at the NMR lab who assisted us with the experiment. I also acknowledge the contribution of Prof. K.R. Shamsundar.

9.7 REFERENCES

- 1. http://www2.chemistry.msu.edu/faculty/reusch/VirtTxtJml/
 Spectrpy/nmr/nmr1.htm
- 2. http://en.wikipedia.org/wiki/Nuclear_magnetic_resonance_ spectroscopy
- 3. http://en.wikipedia.org/wiki/Nuclear_magnetic_resonance
- 4. Fundamentals of NMR, by Thomas L. James
- 5. Long discussion with Prof. K.S. Viswanathan

FLOURESCENCE

November 1, 2012

10.1 OBJECTIVE

To obtain a fluorescence spectrum of the fluorescence of the following kind


- 1. Constant incident frequency, varying detection frequency vs. Intensity detected
- 2. Constant detection frequency (maxima), varying incident frequency vs. Intensity detected
- 3. Varying concentration of solution, with both incident and detection frequencies fixed vs. Intensity detected

10.2 THEORY

Absorption works well when the concentration of the substance to be measured is high. When concentration gets smaller, the absorption gets noise prone and ceases to be reliable. With absorption we can get ppm, but with fluorescence, ppb is a routine affair. However it works IF the molecule fluoresces.

Let's take a look at the fluorescence experiment. We start with a source of light which passes through a diffraction grating, excitation monochromater. Using this, a particular frequency of light is passed through the sample. Now perpendicular to it, there a detector, which comprises of a grating, viz. emission monochromater, and a photomultiplier.

Now before we continue, lets take a deeper look at what's happening. When an incident beam of light strikes a molecule, it goes from the ground electronic state to an excited electronic state. Now depending upon the frequency of the incident beam, the molecule will get into a vibrational state, while it's in the excited electronic state. The vibrational state quickly (within pico seconds typically) drops to the lowest vibrational state, since the molecule is in the liquid state. Consequently, regardless of the incident frequency (granted it's greater than the minimum required to excite the molecule to the first excited electronic state and less then the following electronic state), the molecule will drop from the lowest vibrational energy level in the excited electronic state, to the ground electronic state, at some vibrational state. Thus, if we obtain a graph between the frequency of light

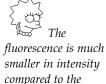

incident beam, thus

perpendicular

	CONCENTRATION (NANO MOLAR)	CPS	
10.0		11070	0.2675
5.0		39355	0.302
2.5		14884	0.635

Table 24: Concentration vs. Intensity

received and intensity, for a fixed incident frequency of light, we'll get a maxima corresponding to the ground vibrational state. Now if use the frequency for which we got the maximum intensity as constant, and change the incident beam's frequency, then again we'll get a maxima. If we use both these maxima values for the incident and emitted beam respectively, and change the concentration, we can use the obtained graph to find unknown concentrations.

We must here realise that there's a phenomenon called Quenching that can kill the emission in fluorescence. Thus, we must perform a time resolved experiment, if a sample were to be analysed for further use.

10.3 OBSERVATIONS/DISCUSSION

The spectrum was not made available. The emission maxima was obtained at 512 nm, and the excitation maxima was obtained at 491 nm. Values for Concentration vs. Intenisty are given in Table 24. The graph is given in Figure 8

10.4 ACKNOWLEDGEMENTS

I thank our PhD assistant who performed the experiment for us, and of course Prof. KS Viswanathan for the rest.

10.5 REFERENCES

1. Prof. K S Viswanathan's Lectures

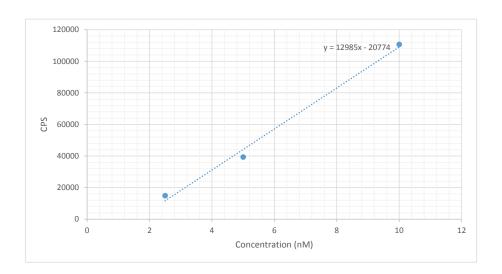


Figure 8

X-RAY DIFFRACTION AND POLARIMETRY

November 8, 2012

11.1 OBJECTIVE

To

- 1. observe and understand the X-ray diffraction technique
- 2. use optical activity of a substance for concentration estimation

11.2 THEORY

11.2.1 X-ray diffraction

The very fist condition for X-ray diffraction analysis of a substance is that it must crystallize. To check the quality of the crystal, a polarized beam of light is passed through it and using a microscope and an analyser (another polariser), the crystal is viewed. The colour must be the same for any given orientation of the analyser, if the crystal formed is good in terms of being analysed by this technique. For holding the crystal, nylon thread loops are used, with super-glue to fix them in place. This is then placed into the single crystal x-ray diffraction machine. There are various degrees of freedom in this device. κ can range from -90° to $+90^{\circ}$ using an arm. The base can rotate 360°. The detector can move 2 θ , where θ depends on the detector being used. Also, the detector can move forward and backwards. Intenisty falls rapidly with distance. However, if the diffracted beams are at a very small angular difference, then distance helps resolve these. An optimal condition must be met to improve resolution at the cost of intensity. The anode is required to be cooled constantly to create the X-ray. To produce higher intensities of X-rays, moving anodes are used, that helps dissipate the heat and thereby allows for higher power. The spots that are obtained, are first matched with a unit cell that's specified by the user. The unit cells can be mathematically shown to be limited. X-ray diffraction is a very strong tool as this is amongst the few analytic techniques that can be used to resolve the chirality of the molecule, which is essentially known if you can resolve the exact 3D-structure of the molecule, which is what is obtained after calculations using this technique. Once the spatial arrangement is known, the R and S classification can be done easily. The temperature is normally kept to the room temperature, but for more

accurate results, the temperature can be dropped to avoid thermal energy induced errors. The machine also has a CO₂ laser which is used to analyse substances that crystallize at low temperatures only. The experimental setup is such that the laser is reflected and targeted using a combination of mirrors, atleast one of which is movable. Now the substance to be analysed is put in a thin high quality capillary tube which is sealed using a super-glue. Now the capillary is placed in the machine and cooled. The laser heats a small part of the wall of the capillary, which in turn heats the now frozen liquid inside back into a liquid. The laser moves upwards gradually, as slow as 3cm in 2 hours. This is similar to zone melting and it results in formation of a crystal. After the process is over, a single transparent frozen area is usually obtained, which is used for analysis. The rest of the area is also solid, but not transparent.

11.2.2 Polarimetry

Specific rotation is given by $\alpha = \theta/lc$, where θ is the angle by which the light is rotated, l is the path length, and c is the concentration given in moles per Litre. Optically active substances are those that rotate the plane of polarization of plane-polarized light. Before we get into understanding what that means, it's important to realise that refractive index of a medium, which essentially dictates the speed of light through the medium, is a function of the material's dipole properties. Having made that clear, let us now decompose our plane polarized light into two components, each consisting of a circularly polarised light, with one going clockwise and the other going anticlockwise, suitably chosen to produce the plane polarized light, when added vectorially. Now if the molecule's environment is chiral, viz. the molecule is not a super-imposable mirror image of itself, on itself¹, and in a given mixture, there exists only one kind of molecule (say right handed using some arbitrary definition), then one of the circularly polarised light (say clockwise) will travel at a different rate than the anti-clockwise circularly polarised light and thereby when it leaves the medium, their superposition would result in a plane wave with a plane that's at an angle with the incident plane of the wave. An important question to address at this stage is that we're talking about molecules which are randomly oriented and what effect will this have on the discussion so far. A little thought reveals that the molecules with the chiral axis at some angle ϕ wrt the plane of light, are equally likely to be found as the molecules at an angle $-\phi$. Since the molecules are chiral, their contribution to the phase difference of the circularly polarised decomposition of the plane polarised light, will not average out to zero. It's the cumulative effect that we see.

 $[\]ensuremath{\mathtt{1}}$ which more precisely may be stated as, only if it does not posess an axis of improper rotation, S_n

From this it's natural to consequence that the optical activity depends directly on the concentration of the molecule. Also, the phase difference will be larger if the light travels a longer distance inside the chiral medium, and thus optical activity depends on the

What we're doing in the experiment is essentially keeping the length constant and changing the concentration. If we plot the θ against the concentration, we will get $l\alpha$, which is constant for the given experiment.

11.3 OBSERVATIONS/DISCUSSION

11.3.1 X-ray diffraction

There aren't any 'observations' to report for this experiment as this experiment was demonstrated to us.

11.3.2 Polarimetry

Values and the graph for Relative Concentration vs. θ are given in Figure 9.

11.4 ACKNOWLEDGEMENTS

I thank our PhD assistant who performed the experiment for us, and of course Prof. KS Viswanathan for the rest.

11.5 REFERENCES

1. Prof. K S Viswanathan's Lectures

Relative Concetration	Optical Rotation	Average Optical Rotation
	0.699	
	0.701	
1	0.7	0.7
	0.281	
	0.282	
	0.283	
0.5	0.283	0.28225
	0.163	
	0.164	
0.25	0.163	0.1633333333

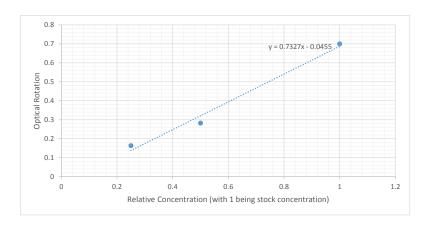


Figure 9

Part II

THE SHOWCASE

You can put some informational part preamble text here. Illo principalmente su nos. Non message *occidental* angloromanic da. Debitas effortio simplificate sia se, auxiliar summarios da que, se avantiate publicationes via. Pan in terra summarios, capital interlingua se que. Al via multo esser specimen, campo responder que da. Le usate medical addresses pro, europa origine sanctificate nos se.

COLOPHON This document was typeset using the typographical look-and-feel classicthesis developed by André Miede, for LATEX. The style was inspired by Robert Bringhurst's seminal book on typography "The Elements of Typographic Style". The latest version of this document is available online at: https://github.com/toAtulArora/IISER_repo